Sensores SpO2

Descripción

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
Instructions for use of Medical Cables SpO2 sensors

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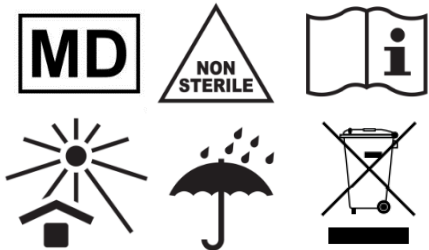medicalcables
Innovation & Compatibility for Health

SpO2 sensors

Description

SpO2 sensors are accessories to be used in conjunction with pulse oximeters or devices for measuring the percentage of oxygen in the blood.

Product performance characteristics – Compatibility

The sensors must be used exclusively with the monitoring equipment indicated on the label on the packaging or on the product label affixed to the sensor cable itself. If the sensor is not compatible with the monitor it is to be used with, the monitor will sound an alarm and display on the screen that the sensor is not working or is not recognized by the device.

If in doubt about product compatibility, consult the manufacturer of each individual device to verify the compatibility of specific models of monitors and sensors.

Our sensors are compatible only in terms of arterial oxygen saturation and pulse rate measurement.

Indications/Contraindications of Medical Cables SpO2 sensors

SpO2 sensors series AC, PC, NW, AE, AS and PS (AC – Adult finger clip; PC – Pediatric finger clip; NW – Neonate silicone wrap; AE – Adult ear clip; AS – Adult soft tip; PS – Pediatric soft tip)

Medical Cables SpO2 sensors series AC, AE and AS are indicated for non-invasive, short-term and continuous monitoring of arterial oxygen saturation and pulse rate in patients over 40 kg body weight. The sensor is contraindicated in patients weighing less than 40 kg and in active patients.

Medical Cables SpO2 sensors series PC and PS are indicated for non-invasive, short-term and continuous monitoring of arterial oxygen saturation and pulse rate in patients less than 40 kg body weight. The sensor is contraindicated in patients weighing more than 40 kg and in active patients.

Medical Cables SpO2 sensors series NW are indicated for non-invasive, short-term and continuous monitoring of arterial oxygen saturation and pulse rate in neonate patients less than 3 kg body weight. The sensor is contraindicated in patients weighing more than 3 kg and in active patients.

SpO2 sensors are not designed for prolonged use or long-term monitoring in one position; it should be repositioned or removed every 2-4 hours, or more frequently if the circulatory status and/or skin integrity so recommend and repositioned at a different site. The sensor must not be fixed at the site of tissue damage.

Instructions for use of Medical Cables SpO2 sensors series AC and PC

- While holding the clip open, place the finger on the bottom pad so that the cable runs along the back of the hand (image 1).
- The tip of the finger should rest against the finger stopper located at the front of the lower pad. The nail should be directly under the light source from the LED (image 2).
- If the chosen finger cannot be positioned correctly, a smaller finger can be used. Do not use the AC and PC series on thumbs or toes, or on a child's hand or foot (image 3).
- When choosing a patient placement site, avoid areas near an arterial catheter, blood pressure cuff, or intravascular infusion line.
- Plug the sensor into the monitoring device and verify the correct operation, according to the instruction manual of the device.
- Before applying the sensor to the patient, the function of the sensor must be tested by trained personnel. When applied to the patient, it must be controlled by specialized personnel.

Instructions for use of Medical Cables SpO2 sensors series AE

- While holding the clip open, place the sensor on patient's earlobe (image 4).
- The LED light source and photoreceptor should make full and complete contact with the skin. Avoid areas that have been pierced.
- Do not use SpO2 sensors series AE on any other site than the patient ears.
- Plug the sensor into the monitoring device and verify the correct operation, according to the instruction manual of the device.
- Before applying the sensor to the patient, the function of the sensor must be tested by trained personnel. When applied to the patient, it must be controlled by specialized personnel.

Instructions for use of Medical Cables SpO2 sensors series AS and PS

- While holding the soft tip open and so that the finger symbol is at the top, place your finger inside the finger cot, with the cable passing through the back of the hand (images 5 and 6).
- The bottom of the fingernail should be positioned under the LED light.
- If the chosen finger cannot be positioned correctly, it can be placed on a smaller finger. Do not use the AS and PS series on thumbs or toes, or on a child's hand or foot.
- When choosing a patient placement site, avoid areas near an arterial catheter, blood pressure cuff, or intravascular infusion line.
- Plug the sensor into the monitoring device and verify the correct operation, according to the instruction manual of the device.
- Before applying the sensor to the patient, the function of the sensor must be tested by trained personnel. When applied to the patient, it must be controlled by specialized personnel.

Instructions for use of Medical Cables SpO2 sensors series NW

- While holding the neonate wrap open, place the sensor around the patient's foot (image 7).
- The LED light source and photoreceptor should make full and complete contact with the skin.
- Do not use SpO2 sensors series NW on any other site than the patient foots.
- Plug the sensor into the monitoring device and verify the correct operation, according to the instruction manual of the device.
- Before applying the sensor to the patient, the function of the sensor must be tested by trained personnel. When applied to the patient, it must be controlled by specialized personnel.

NOTE: If the monitoring device does not show a reliable reading, it may be because the sensor is not positioned correctly or the sensor site is too thick, thin, or deeply pigmented to allow light to pass through properly. Reposition the sensor or move it to another site.

Medical Cables SpO2 sensors should always be used by healthcare professionals caring for patients with good or poor perfusion, in hospitals and hospital-type facilities.

Clinical benefits of Medical Cables SpO2 sensors

- Continuous oxygenation monitoring: facilitates continuous monitoring of the patient's blood oxygen saturation, allowing a constant evaluation of their oxygenation status through the monitoring equipment to which it is connected.
- Early detection of hypoxemia: contributes to the early detection of low blood oxygen levels by transmitting precise signals to the monitoring equipment, alerting healthcare personnel to potential respiratory or cardiovascular problems.
- Therapy optimization: enables more precise therapy by providing real-time signals to the monitor, helping to adjust the delivery of supplemental oxygen or other therapeutic interventions as needed.
- Treatment response monitoring: facilitates monitoring of the patient's response to treatment by transmitting continuous signals to the monitoring equipment, allowing the effectiveness of therapeutic interventions to be evaluated over time.
- Monitoring during medical procedures: helps lung function monitoring during medical procedures, such as mechanical ventilation, by sending precise signals to the monitoring equipment, helping to optimize therapy and prevent respiratory complications.

Cleaning and disinfection of all models of Medical Cables SpO2 sensors

SpO2 sensors must be cleaned and disinfected after each use, before use on another patient.

To clean and disinfect the sensors:

- Disconnect the sensor from the monitoring device.

- Wet a clean, dry cloth or gauze pad with a mild, low-level disinfectant solution, such as a 1:10 bleach solution. Do not use undiluted bleach as this may cause permanent damage to the sensor. It is used to clean all exposed surfaces of the sensor and cable. Follow the disinfectant manufacturer's instructions for use.
- Wet another clean, dry cloth or gauze pad with distilled water and re-clean all exposed surfaces of the sensor and cable.
- Dry the sensor and cable with a clean dry cloth or gauze.

Medical Cables SpO2 sensors are reusable and are not limited in number of uses; they can always be reused, until they break.

Caution: Do not immerse the sensor in any liquid or use abrasive or corrosive liquids. Do not sterilize by irradiation, steam or ethylene oxide (ETO). Not suitable for autoclaving.

WARNINGS

- All sensors are designed for use with specific monitors. Before use, check the compatibility of the monitor and the sensor, otherwise this may cause device performance failure and/or harm to the patient.
- The sensor shall be free from visible defects, discoloration and damage. If the sensor is discolored or damaged, discontinue use. Never use a damaged sensor or one that has damaged electrical circuits.
- If the SpO2 sensor is not positioned properly it may give incorrect measurements.
- Nail varnish or artificial nails can cause inaccurate readings, so they must be removed before placing the sensor. The same can happen in patients with dark or highly pigmented skin.
- Reusable sensors should be repositioned periodically as mentioned above (see indications/contraindications). Because skin tolerance to the sensor varies from patient to patient, some people require more frequent sensor site changes
- Intravascular contrast media can lead to inaccurate measurements.
- Movement can affect the function of the SpO2 sensors (therefore use on active patients is contraindicated). It is also contraindicated in patients with low perfusion.
- Do not use other than dedicated tapes to hold the sensor in place as the additional pressure may cause venous pulsations which can lead to inaccurate measurements.
- Be sure to position the sensor cable correctly to reduce the chance of it becoming entangled in the patient, which could damage the sensor, cable, or monitoring device.
- Intense sources of visible and infrared light, such as direct sunlight, bilirubin lights, or infrared heaters can interfere with the proper function of the SpO2 sensors and cause inaccurate readings. If this situation occurs, cover the sensor with an opaque cloth.
- Proper function of the SpO2 sensor may be affected by electromagnetic sources such as electrosurgical equipment.
- Do not use the sensor during Magnetic Resonance Imaging (MRI) or Computed Axial Tomography (CAT) procedures. Conducted current can cause burns and the MR or CT unit can cause inaccurate readings; and the sensor can affect the MR or CT images.
- Do not pull on the sensor cable when disconnecting it from the device, as this may damage the sensor and device. Always hold it by the connector and remove it carefully.
- If the cable or sensor is damaged in any way, stop using it immediately and replace it with another sensor of the same type.
- Do not alter or modify the sensor in any way as this may affect its performance and accuracy and will void the product warranty.
- This SpO2 sensor is for use in protection class II/Type BF or CF medical equipment only.

Warranty

Products come with our manufacturer’s warranty.

Medical Cables SpO2 sensors are warranted against defects in material and/or workmanship for a period of 3 years from the date of shipment. During this warranty period, Medical Cables will replace, free of charge, any product found to be defective.

This warranty will not be valid in the event of misuse, improper use of the product, improper handling of the product by an unauthorized person, deterioration of the product due to external agents such as water or other harmful chemical or corrosive products, inadequate maintenance or lack of cleaning, or any fault that is the result of normal wear and tear. Alterations or modifications to the product void this warranty.

The conformity of the products is guaranteed according to the use for which they are intended.

Return policy

The term for changes or returns of purchased products is 12 months from the date of shipment of the order. It is required to notify your sales representative by email mentioning the reason for the change or return. The commercial will assess your request and will contact you to arrange the collection of the material.

Technical data

Measuring range: 0-100%
Pulse rate range: 20-250 BPM
Resolution: SpO2 ±1%; Heart rate ±1BPM
Accuracy: SpO2 70%-100% ±3%; unspecified <69%. Heart rate ±2BPM.

Specific storage and handling conditions

Ideally, products should be stored in the same wrapper or packaging in which they were supplied by the manufacturer.

Operating temperature: 5°C ~ 40°C
Operating relative humidity: 30-75%
Storage temperature: -20°C ~ 55°C
Storage humidity: <80% non-condensing

Waste

Do not throw garbage. Check local laws and regulations for information on how to dispose of SpO2 sensors.

WARNING

Any serious incident related to the product must be reported to the manufacturer and to the competent authority of the Member State in which the user and/or the patient are established.

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